

LIFEPAK[®] 15

monitor/defibrillator



For emergency medical services



When you respond to emergencies, you need the most advanced monitor/defibrillator that sets the standard in innovation, operations and toughness.



The LIFEPAK 15 monitor/defibrillator delivers


Physio-Control defibrillators have set the standard for over 60 years, and the enhanced LIFEPAK 15 monitor/defibrillator raises the bar. As our most advanced emergency response monitor/defibrillator, the LIFEPAK 15 device balances sophisticated clinical technologies and supreme ease of use in a device that's tough enough to stand up to your most challenging environments. Evolving from its original platform, the LIFEPAK 15 adds new features—temperature monitoring and external power—to complement existing features which include 360J energy and 12-lead ECG transmission. And that means your team can be even more effective.

A LIFEPAK device never stands on its own—and the LIFEPAK 15 monitor is no different. Physio-Control is committed to providing innovative solutions for emergency response care, from first responders to throughout the hospital.

Our products have helped save tens of thousands of lives. We're proud to continue this work with new features in the LIFEPAK 15 monitor/defibrillator.

The standard in clinical innovation

The pioneer in portable defibrillation and monitoring technology, Physio-Control is committed to creating technologies and devices that change the way you provide emergency care. You can see the results in the enhanced LIFEPAK 15 monitor/defibrillator, which sets a new standard in innovation—yet again.

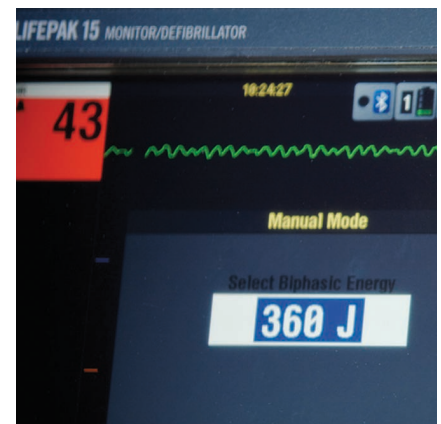

Advanced monitoring parameters

The LIFEPAK 15 gives you more monitoring capabilities than any other device. Noninvasive monitoring of carbon monoxide, SpO₂, and methemoglobin (related to certain chemical exposures and drugs) offered



by Masimo® Rainbow® Technology helps you detect hard-to-diagnose conditions and improve patient care. In addition, the 15 now offers temperature monitoring—and like other data, you can transmit it to other systems, trend it, or display for post-event review in CODE-STAT™ software.

Advanced support for treating cardiac patients

The LIFEPAK 15 continuously monitors all 12-leads in the background and alerts you to changes using the ST-Segment trend monitoring feature, after acquiring the initial 12-lead. Additionally, STJ values are now included on the 12-lead printout to help you identify changes. The LIFEPAK 15 also works seamlessly with the web-based LIFENET System 5.0, so you can automatically share critical patient data with multiple patient care teams.

Full energy up to 360 joules, for every patient who needs it

The LIFEPAK 15 monitor/defibrillator features 360J biphasic technology, which gives you the option of escalating your energy dose up to 360J for difficult-to-defibrillate patients. Why is this necessary? Recent studies have shown that refrillation is common among VF cardiac arrest patients and that defibrillation of recurring episodes of VF is increasingly difficult. A randomized controlled clinical trial shows the rate of VF termination was higher with an escalating higher energy regimen of 200J and over.¹

Proven CPR guidance and post event review

The CPR Metronome in the LIFEPAK 15 monitor uses audible prompts to guide you without distracting vocal critique. A metronome has been a feature that has been demonstrated to help professionals perform compressions and ventilations within the recommended range of the 2015 ERC Guidelines. Post-event review of CPR data and delivering feedback to the team has been shown to be effective in improving CPR quality in both hospital and out-of-hospital.^{2,3,4} And by transmitting code data directly to CODE-STAT Data Review software, EMS personnel can review CPR statistics and provide training and feedback where it is most needed.



The standard in operational effectiveness

Flexible, connected and easy to use, the LIFEPAK 15 monitor/defibrillator was designed based on the feedback and needs specific to working in the field.

Dual-mode LCD screen with SunVue™ display

Switch from full-color to high-contrast SunVue mode with a single touch for the best full-glare view in the industry. A large screen (8.4 inches diagonally) and full-color display provide maximum viewability from all angles.

Flexible power options

Choose between external worldwide AC or DC power, or use the latest Lithium-ion dual battery technology for up to six hours of power. The LIFEPAK 15 monitor's two-battery system requires no maintenance or conditioning, and allows you to charge batteries in the device. In addition, you can track the status and service life of your batteries using LIFENET® Asset, part of the LIFENET System data network.

Data connectivity

The LIFEPAK 15 collects code summaries and equipment status data along with critical clinical information as you treat patients. Using LIFENET Connect, part of the LIFENET System data network, the code summaries can be sent directly to your quality improvement team for review with CODE-STAT Data Review Software. Your equipment manager can also view equipment status on the LIFENET System using LIFENET Asset and alert you to any potential issues.

Upgradable platform

The LIFEPAK 15 platform is flexible enough to adapt to evolving protocols and new guidelines, and can be upgraded as you're ready to deliver new capabilities. With more processing power and speed, the LIFEPAK 15 is designed to grow as your needs change, helping you avoid costly premature replacements.

Attention to detail

The LIFEPAK 15 monitor is designed based on field feedback to make it a more effective tool. The LIFEPAK 15 has a larger handle for easier handoffs, an easy to clean keypad, and a common interface to the LIFEPAK 12 defibrillator/monitor that helps reduce training.



LIFEPAK 15 MONITOR/DEFIBRILLATOR

Recommended Adult VF Dose: 200-300 J

PHYSIO CONTROL



SpO2
Calmson SET
SpO2

12-LEAD
TRANSMIT
PRINT

ON
ENERGY SELECT
CHARGE
PACER
RATE
CURRENT
USE

in the presence of flammable gases.
For use only by qualified personnel.

The standard in toughness

We believe LIFEPAK equipment should live up to the highest expectations of those working in the harshest settings. The LIFEPAK 15 is LIFEPAK TOUGH™, with improved ruggedness and durability you can rely on.

Works when dropped, kicked, soaked or dirty

The LIFEPAK 15 monitor/defibrillator passes 30-inch drop tests, which is equal to falling off a cot or dropping it in transit. And with an IP44 rating, it doesn't matter how wet or dirty it gets, so you can keep working in steady wind, rain and other harsh environments.

Toughened inside and out

We heard from emergency response teams that they wanted a tougher device—so we added a shock-absorbing handle, a double-layer screen that can take a beating from doorknobs and cot handles, and redesigned cable connections for confident monitoring and therapy delivery.

Unmatched field service

The unit's self-checking feature alerts our service team if the device needs attention. Our onsite maintenance and repair, access to original manufacturer parts, and highly trained, experienced service representatives give you the peace of mind that your LIFEPAK 15 monitor will be ready when you need it.*



Data connectivity



LIFEPAK TOUGH™



Dual-mode LCD screen with SunVue display

*A variety of customized service options are available.




LIFEPAK 15 monitor/defibrillator

The latest lithium-ion battery technology and dual battery system allows for nearly six hour run time, automatic switching between external power and batteries, and an approximate two-year replacement cycle.

Easy one-touch Bluetooth® data transmission.

12-lead ECG transmissions via the LIFENET System and ST-Segment trend monitoring make the LIFEPAK 15 unit a vital part of decreasing EMS-to-balloon (E2B) response times.

Integrated carbon monoxide and methemoglobin monitoring.

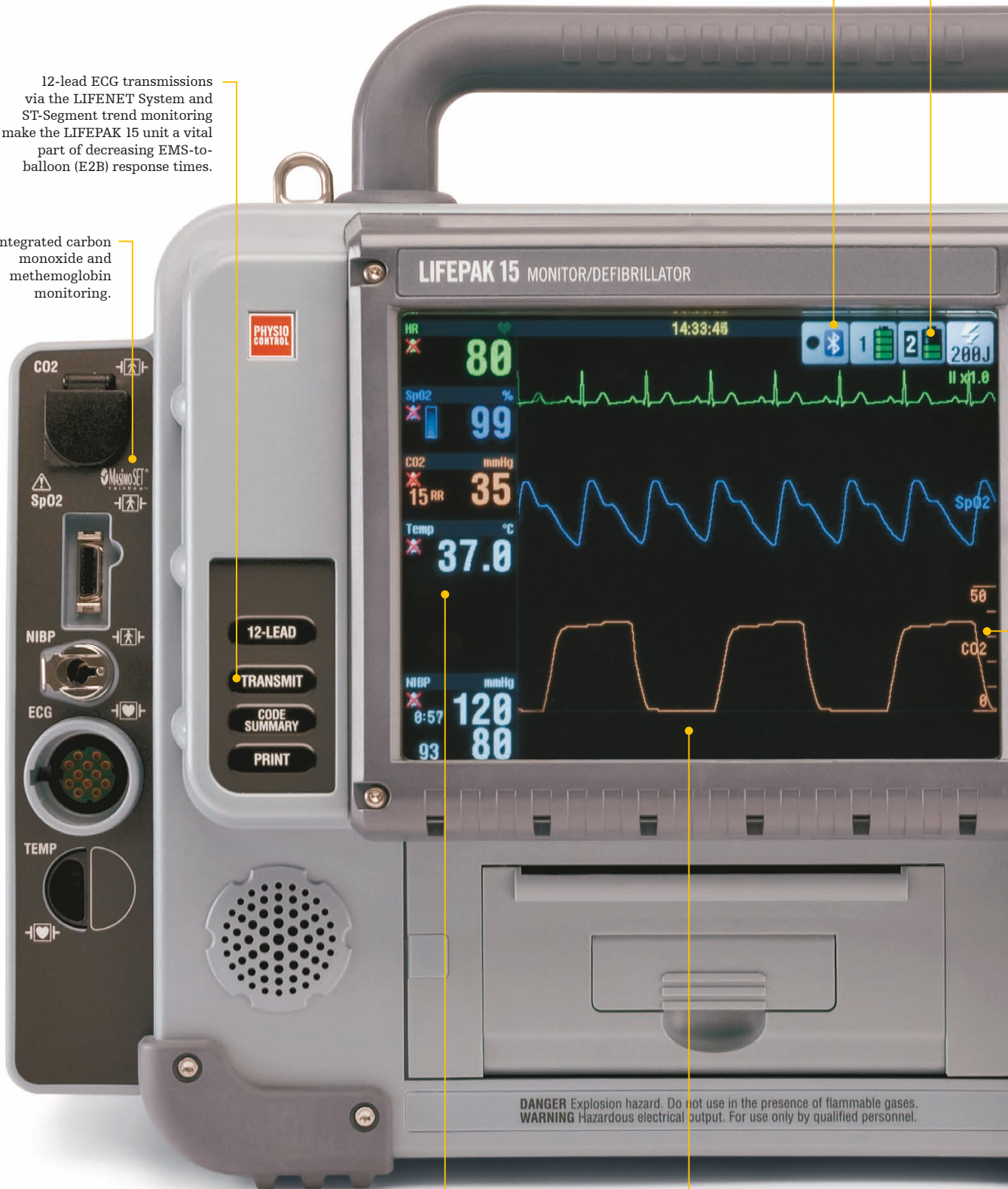

On-screen temperature display in either Celsius or Fahrenheit.

Large screen with high resolution for maximum viewability in any setting.

Ergonomically designed handle has built-in shock absorbers for cushion and fits two gloved hands for easy pass off.



CPR Metronome is a proven technology that actively guides users to a consistent compression rate without the need for extra external hardware.

Integrated Microstream™ EtCO₂ monitoring by Medtronic provides waveform ranges as low as 0-20 mmHg to help identify ROSC or gauge CPR quality, consistent with current ERC guidelines.

The LIFEPAK 15 monitor/ defibrillator at a glance

Redesigned cable connector gives you the confidence for secure therapy delivery.



For more than 60 years, Physio-Control has been developing technologies and designing devices that are legendary among first response professionals, clinical care providers, and the community.



A legacy of trust

Since we were founded in 1955, Physio-Control has been giving medical professionals around the world legendary quality and constant innovation. Our LIFEPAK devices have been carried to the top of Mount Everest. They've been launched into orbit on the International Space Station. And you'll find more than half a million units in use today on fire rescue rigs, ambulances, and hospital crash carts worldwide.

We are inspired and informed by the rescuers who choose our products to save lives. The knowledge gained from working with some of the world's largest EMS organizations helps us constantly improve clinical standards and durability.

Today, we continue our legacy of innovation with leading technologies that improve patient care. Our 360J biphasic technology gives patients the best chance at survival. Our secure, web-based flow of ECG data helps improve STEMI patient outcomes. And our carbon monoxide monitoring helps catch the number one cause of poisoning deaths.

From the streets to the emergency room to the administrative office, we offer a powerful suite of solutions that range from code response to quality control analysis. And even as we bring ground-breaking products to the market, some things don't change. As always, when you choose our products, you don't just get a device. You also get the most comprehensive warranty in the business, industry-leading technical service, and a partner with over 60 years of experience in emergency care.

For more information about the enhanced LIFEPAK 15 monitor/defibrillator—and how it can help you do what you do best—please contact your local Physio-Control representative or visit www.physio-control.com.

Defibrillator

Biphasic waveform: Biphasic Truncated Exponential

The following specifications apply from 25 to 200 ohms, unless otherwise specified:

Energy accuracy: ± 1 joule or 10% of setting, whichever is greater, into 50 ohms, ± 2 joules or 15% of setting, whichever is greater, into 25-175 ohms.

Voltage compensation: Active when disposable therapy electrodes are attached. Energy output within $\pm 5\%$ or ± 1 joule, whichever is greater, of 50 ohms value, limited to the available energy which results in the delivery of 360 joules into 50 ohms.

Paddle options: QUIK-COMBO pacing/defibrillation/ECG electrodes (standard). Cable Length 8 foot long (2.4 m) QUIK-COMBO cable (not including electrode assembly). Standard paddles (optional)

Manual mode

Energy select: 2, 3, 4, 5, 6, 7, 8, 9, 10, 15, 20, 30, 50, 70, 100, 125, 150, 175, 200, 225, 250, 275, 300, 325, and 360 joules

Charge time: Charge time to 360 joules in less than 10 seconds, typical

Synchronous cardioversion: Energy transfer begins within 60 msec of the QRS peak

Paddles lead off sensing: When using QUIK-COMBO electrodes, the device indicates Paddles Leads Off if the resistive part of the patient impedance is greater than $300 \pm 15\%$ ohms, or if the magnitude of the patient impedance is greater than $440 \pm 15\%$ ohms.

AED mode

Shock Advisory System™ (SAS): an ECG analysis system that advises the operator if the algorithm detects a shockable or non-shockable ECG rhythm. SAS acquires ECG via therapy electrodes only.

Shock ready time: Using a fully charged battery at normal room temperature, the device is ready to shock within 20 seconds if the initial rhythm finding is "SHOCK ADVISED"

Biphasic output: Energy Shock levels ranging from 150–360 joules with same or greater energy level for each successive shock

cprMAX™ Technology: In AED mode, cprMAX technology provides a method of maximizing the CPR time that a patient receives, with the overall goal of improving the rate of survival of patients treated with AEDs.

Setup options:

- Auto Analyze: Allows for auto analysis. Options are OFF, AFTER 1ST SHOCK
- Initial CPR: Allows the user to be prompted for CPR for a period of time prior to other activity. Options are OFF, ANALYZE FIRST, CPR FIRST
- Initial CPR Time: Time interval for Initial CPR. Options are 15, 30, 45, 60, 90, 120, and 180 seconds.
- Pre-Shock CPR: Allows the user to be prompted for CPR while the device is charging. Options are OFF, 15, 30 seconds.
- Pulse Check: Allows the user to be prompted for a pulse check at various times. Options are ALWAYS, AFTER EVERY SECOND NSA, AFTER EVERY NSA, NEVER
- Stacked Shocks: Allows for CPR after 3 consecutive shocks or after a single shock. Options are OFF, ON
- CPR Time: 1 or 2 User selectable times for CPR. Options are 15, 30, 45, 60, 90, 120, 180 seconds and 30 minutes.

Pacer

Pacing mode: Demand or non-demand rate and current defaults

Pacing rate: 40 to 170 PPM

Rate accuracy: $\pm 1.5\%$ over entire range

Output waveform: Monophasic, truncated exponential current pulse (20 ± 1 ms)

Output current: 0 to 200 mA

Pause: Pacing pulse frequency reduced by a factor of 4 when activated

Refractory period: 180 to 280 msec (function of rate)

Environmental

Unit meets functional requirements during exposure to the following environments unless otherwise stated.

Operating temperature: 0° to 45°C (32° to 113°F); -20°C (-4°F) for 1 hour after storage at room temperature; 60°C (140°F) for 1 hour after storage at room temperature

Storage temperature: -20° to 65°C (-4° to 149°F) except therapy electrodes and batteries

Relative humidity, operating: 5 to 95%, non-condensing. NIBP: 15 to 95%, non-condensing

Relative humidity, storage: 10 to 95%, non-condensing

Atmospheric pressure, operating: -382 to 4,572 m (-1,253 to 15,000 ft). NIBP: -152 to 3,048 m (-500 to 10,000 ft)

Water resistance, operating: IP44 (splash proof, dust and sand resistant) per IEC 529 and EN 1789 (without accessories except for 12-lead ECG cable, hard paddles, and battery pack)

Vibration: MIL-STD-810E Method 514.4, Propeller Aircraft - category 4 (figure 514.4-7 spectrum a), Helicopter - category 6 (3.75 Grms), Ground Mobile - category 8 (3.14 Grms), EN 1789: Sinusoidal Sweep, 1 octave/min, 10-150 Hz, ± 0.15 mm/2 g

Shock (drop): 5 drops on each side from 18 inches onto a steel surface EN 1789: 30-inch drop onto each of 6 surfaces

Shock (functional): Meets IEC 60068-2-27 and MIL-STD-810E shock requirements 3 shocks per face at 40 g, 6 ms half-sine pulses

Bump: 1000 bumps at 15 g with pulse duration of 6 msec

Impact, non-operating: EN 60601-1 0.5 + 0.05 joule impact UL 60601-1 6.78 Nm impact with 2-inch diameter steel ball. Meets IEC62262 protection level IK 04.

EMC: EN 60601-1-2:2006 Medical Equipment - General Requirements for Safety - Collateral Standard: Electromagnetic Compatibility - Requirements and Tests EN 60601-2-4:2003: (Clause 36) Particular Requirements for the Safety of Cardiac Defibrillators and Cardiac Defibrillator-Monitors

Cleaning: Cleaning 20 times with the following: Quaternary ammonium, isopropyl alcohol, hydrogen peroxide

Chemical resistance: 60 hour exposure to specified chemicals: Betadine (10% Povidone-Iodine solution), Coffee, Cola, Dextrose (5% Glucose solution), Electrode Gel/Paste (98% water, 2% Carbopol 940), HCL (0.5% solution, pH=1), Isopropyl Alcohol, NaCl solution (0.9% solution), Cosmetic discoloration of the paddle well shorting bar shall be allowed following exposure to HCL (0.5% solution).

Power

Power adapters: AC or DC

Power Adapters provide operation and battery charging from external AC or DC power

- Full functionality with or without batteries when connected to external AC/DC
- Typical battery charge time while installed in LIFEPAK 15 device is 190 minutes
- Indicators: external power indicator, battery charging indicator

Dual battery: Capability with automatic switching

Low battery indication and message: Low battery fuel gauge indication and low battery message in status area for each battery

Replace battery indication and message: Replace battery fuel gauge indication, audio tones and replace battery message in the status area for each battery. When replace battery is indicated, device auto-switches to second battery. When both batteries reach replace battery condition, a voice prompt instructs user to replace battery.

Battery capacity

For two, new fully-charged batteries, 20°C (68°F)

Operating mode	Monitoring (minutes)	Pacing (minutes)	Defibrillation (360J discharges)	
			Typical	Minimum
Total capacity to shutdown	Typical	360	340	420
	Minimum	340	320	400
Capacity after low battery	Typical	21	20	30
	Minimum	12	10	6

Battery

Battery specifications

Battery type: Lithium-ion

Weight: < 0.60 kg (1.3 lb)

Charge time (with fully depleted battery): 4 hours and 15 minutes (typical)

Battery indicators: Each battery has a fuel gauge that indicates its approximate charge. A fuel gauge indicates that the battery should be replaced.

Charging temperature range: 5° to 45°C (41° to 113°F)

Operating temperature range: 0° to 45°C (32° to 113°F)

Short term (<1 week) storage temperature

Range: -20° to 60°C (-4° to 140°F)

Long term (>1 week) storage temperature

Range: 20° to 25°C (68° to 77°F)

Operating and storage humidity range: 5 to 95% relative humidity, non-condensing

References

- 1 Stiel I, Walker R, Nesbitt L, et al. Biphasic Trial: A randomized comparison of fixed lower versus escalating higher energy levels for defibrillation in out-of-hospital cardiac arrest. *Circulation*. 2007;115:1511-1517.
- 2 Edelson D, Litzinger B, Arora V, et al. Improving in-hospital cardiac arrest process and outcomes with performance debriefing. *Arch Intern Med*. 2008;168:1063-1069.
- 3 Olasveengen T, Wik L, Kramer-Johansen J, et al. Is CPR quality improving? A retrospective study of out-of-hospital cardiac arrest. *Resuscitation*. 2007;75:260-266.
- 4 Fletcher D, Galloway R, Chamberlain D, et al. Basics in advanced life support: A role for download audit and metronome. *Resuscitation*. 2008;78:127-134.

Physio-Control Headquarters
11811 Willows Road NE
Redmond, WA 98052 USA
Tel 425 867 4000
www.physio-control.com

 **Physio-Control, Inc.**, 11811 Willows Road NE, Redmond, WA 98052 USA

 0123

Siège: 27, Old Gloucester Street, LONDON, WC1N 3AX, (UK) - Tél. +44 20 3398 1339
Bureaux: 15 rue de la Garenne 27950 Saint-Marcel (France) Tél. +33 1 87 66 35 60
<https://ibsgroupe.com> – Email: dircom@ibsgroupe.com

Turnkey Supplier